

This discussion paper is/has been under review for the journal Atmospheric Measurement Techniques (AMT). Please refer to the corresponding final paper in AMT if available.

# Quantitative and enantioselective analysis of monoterpenes from plant chambers and in ambient air using SPME

**N. Yassaa<sup>1,2</sup>, T. Custer<sup>1</sup>, W. Song<sup>1</sup>, F. Pech<sup>3</sup>, J. Kesselmeier<sup>3</sup>, and J. Williams<sup>1</sup>**

<sup>1</sup>Air Chemistry Department, Max-Planck Institute for Chemistry, J.J. Becher Weg 27, 55020 Mainz, Germany

<sup>2</sup>Faculty of Chemistry, University of Sciences and Technology Houari Boumediene, U.S.T.H.B., B.P. 32 El-Alia, Bab-Ezzouar, 16111 Algiers, Algeria

<sup>3</sup>Biogeochemistry Department, Max-Planck Institute for Chemistry, J.J. Becher Weg 27, 55020 Mainz, Germany

Received: 20 July 2010 – Accepted: 27 July 2010 – Published: 9 August 2010

Correspondence to: N. Yassaa (nyassaa@usthb.dz)

Published by Copernicus Publications on behalf of the European Geosciences Union.

**AMTD**

3, 3345–3381, 2010

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

## Abstract

A solid-phase microextraction (HS-SPME) and gas chromatography/mass spectrometry (GC/MS) system has been developed for quantifying enantiomeric and nonenantiomeric monoterpenes in plant chamber studies and ambient air. Performance of this system was checked using a capillary diffusion system to produce monoterpene standards. The adsorption efficiency, competitive adsorption and chromatographic peak resolution of monoterpene enantiomer pairs were compared for three SPME fibre coatings: 75  $\mu\text{m}$  Carboxen-PDMS (CAR-PDMS), 50/30  $\mu\text{m}$ , divinylbenzene-carboxen-polydimethylsiloxane (DVB-CAR-PDMS) and 65  $\mu\text{m}$  divinylbenzene-polydimethylsiloxane (DVB-PDMS). Key parameters such as the linearity and reproducibility of the SPME system have been investigated in this work. The best compromise between the enantiomeric separation of monoterpenes and competitive adsorption of the isoprenoids on the solid SPME fibre coating was found for DVB-PDMS fibres. The optimum conditions using DVB-PDMS fibres were applied to measure the exchange rates of monoterpenes in the emission of *Quercus ilex* using a laboratory whole plant enclosure under light and dark conditions, as well as in ambient air. With 592 and 223  $\text{ng m}^{-2} \text{s}^{-1}$ , respectively,  $\beta$ -myrcene and limonene were the predominant monoterpenes in the emission of *Q. ilex*. These values were closely comparable to those obtained using a zNose and cartridge GC-FID systems.

## 1 Introduction

Biogenic emissions of volatile organic compounds (VOCs) are estimated to exceed anthropogenic emissions by roughly a factor of ten on the global scale (Guenther et al., 1995; Müller, 1992). Due to the high reactivity of BVOC with the primary oxidizing agent OH and their propensity to form particles, these emissions can affect global oxidation chemistry and the Earth's radiation budget (Williams, 2004 and references therein). One much studied sub-set of BVOC is the family of monoterpenes, which are emitted

AMTD

3, 3345–3381, 2010

### Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

**Quantitative analysis  
of chiral  
monoterpenes using  
SPME**

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



by terrestrial plants to the atmosphere in large amounts (ca.  $150 \text{ TgC y}^{-1}$ ,  $T_g = 10^{12} \text{ g}$ ) (Guenther et al., 1995). Monoterpenes are particularly important in atmospheric chemistry, reacting rapidly with ozone,  $\text{NO}_3$  and OH, often leading to the production of secondary organic particles in significant yield. Their function in the natural world may be as a signaling compound for insects, aiding pathogen and pollinators alike in finding plants. They might also act as a chemical defense against various insects as well as play a role in inter plant communication (Croteau, 1987; Baldwin et al., 2006; Unsicker et al., 2009). Most of these monoterpenes are chiral, although in most atmospheric chemistry studies the (+)- and (–)-enantiomers are not resolved and hence measured together (e.g., Kesselmeier and Staudt, 1999). While the gas phase reaction rates of both enantiomers are identical, the individual enantiomers often have markedly different biological properties (Lough, 2002). Interestingly, it has been recently reported that optical isomers can react at significantly different rates when they are adsorbed on the surface of aerosol (Stokes et al., 2009).

Reliable and quantitative methods are needed for the analysis of chiral monoterpenes in order to evaluate to what extent different optical isomers contribute to formation of secondary organic aerosols (SOA), to understand the metabolic processes underlying their release, and to study their ambient distribution. To date, little work has been done to analyze naturally occurring enantiomeric gas phase biogenic emissions. Specific plant emissions have been investigated previously using a  $\beta$ -cyclodextrin capillary GC with charcoal adsorption/solvent extraction (Yassaa et al., 2000), and carbograph adsorption/thermal desorption (Yassaa et al., 2001; Williams et al., 2007). In our previous works, we have reported on the applicability of SPME for qualitative study of the biogenic emission of enantiomeric and non-enantiomeric monoterpenes from living plants (Yassaa and Williams, 2005, 2007). For the first time, we present a quantitative method to measure enantiomeric and nonenantiomeric monoterpenes using SPME to extract monoterpenes. The SPME system was calibrated using standards generated by a home-built capillary diffusion system. Three different SPME fibre coatings were tested, including  $75 \mu\text{m}$  Carboxen-PDMS (CAR-

**Quantitative analysis  
of chiral  
monoterpenes using  
SPME**

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



PDMS), 50/30  $\mu\text{m}$  divinylbenzene-carboxen-polydimethylsiloxane (DVB-CAR-PDMS), and 65  $\mu\text{m}$  divinylbenzene-polydimethylsiloxane (DVB-PDMS). Each of these coatings were compared with regard to their adsorption efficiency, competitive adsorption effects and achievable chromatographic peak resolution for pairs of monoterpene enantiomers. Key parameters such as the evolution of the extraction efficiencies as a function of sampling time, the linearity of SPME adsorption and detection limit were also investigated in the course of this work. An optimised method was then employed to measure enantiomeric and nonenantiomeric monoterpenes emitted by *Quercus ilex* in laboratory whole plant enclosure as well as in ambient air. A comparison of simultaneous measurements of these emissions using SPME, solid sorbent GC-FID and a zNose (fast chromatograph equipped with a surface acoustic wave sensor) provide further validation of this method.

## 2 Experimental

### 2.1 SPME

#### 2.1.1 Standard generation

Monoterpene standards were generated using a 5 channel capillary diffusion system (CDS) as displayed in Fig. 1, which employs similar principle as that used by e.g., Helmig et al. (2003) and Baker and Sinnott (2009). Briefly, this system consists of an aluminium block fitted with 2 peltier elements for temperature regulation (maintained for these experiments at  $27 \pm 0.02^\circ\text{C}$  using custom built PID regulation system) and drilled through to house five glass chambers. Glass chambers (160 mm length, 2.22 mm o.d., 1.9 mm i.d.) are fitted with custom built Teflon end caps accommodating insertion of 0.32 mm (1/8 inch) tubing. A regulated flow ( $\sim 50$  SCCM, see Table 1, maintained by MKS mass flow controllers) of compressed air was passed through two ozone scrubbers (Pollman et al., 2005) to avoid any oxidation of monoterpene species

and then through each of these chambers independently. Approximately 500  $\mu\text{L}$  of selected pure monoterpene standards (Sigma–Aldrich, Germany) were placed separately in vials (5 mm quartz NMR tubes, Kontes glass company 897193-0050) and capped using standard NMR caps. The vials were then placed individually into the larger glass chambers, see Fig. 1. A deactivated diffusion capillary (0.53 mm i.d., 5 cm length, Phenyl-Sil-Deactivation) was fed through the cap of each vial to allow a limited amount of chemical to diffuse into the air flowing through each chamber. Output of the CDS was further diluted by mixing the output of the system with a stream of clean compressed air whose flow rate was controlled by a mass flow controller. The weight loss of chemical over time was monitored every few weeks, thus allowing calculation of the theoretical mixing ratio emanating from the CDS. The output of the CDS was also directly measured and compared to measurements of a trace gas standard (National Physical Laboratory, Teddington Middlesex, UK). Results of this comparison are given in Table 1.

### 2.1.2 SPME sampling chamber

For SPME sampling (see Fig. 1), a small septum (Sigma–Aldrich, Supelco, Germany) as well as two 0.64 cm ( $\frac{1}{4}$  inch) glass ports (inlet and outlet) with teflon stopcocks were affixed to a 5 L glass chamber (Skötk, Mainz, Germany). The 5 L volume of the chamber was sufficient to homogeneously mix the analytes prior to SPME sampling. Mixing was initially tested using two 250 mL gas sampling bulbs (Supelco) individually or in series. When sampling from the CDS or standard gas bottle, upon changing chemical mixing ratios the system was allowed to reach steady state by allowing gas to flow through the system for at least an hour. After a steady state was attained, the stopcocks were closed isolating the collected, equilibrated air. The monoterpenes were sampled from the glass chamber, by piercing the septum and exposing the fibre to the equilibrated air for a measured time. Temperature during measurements remained essentially constant at 25 °C in the climate controlled laboratory. Immediately following exposure, the fibres were analysed by GC-MS. The time between SPME sampling and

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



GC-MS analysis never exceeded 3 min. Each SPME sampling was repeated at least 3 times. After SPME sampling, the system was again flushed by re-opening the two stopcocks.

### 2.1.3 Comparison of SPME fibre coatings

5 Three SPME fibre coatings were selected for testing and were all obtained from Supelco (Taufkirchen, Germany). The fibres tested included 75  $\mu\text{m}$  Carboxen-PDMS (CAR-PDMS), 50/30  $\mu\text{m}$ , divinylbenzene-carboxen-polydimethylsiloxan (DVB-CAR-PDMS) and 65  $\mu\text{m}$  divinylbenzene-polydimethylsiloxane-(DVB-PDMS). These were chosen based on previous experience with these polymers (Yassaa and Williams,  
10 2005, 2007) and because their behavior has been well documented in the past (e.g., Bicchi et al., 2000). Each new fibre was thermally conditioned prior to use according to the manufacturer's recommendations. For this comparison, each of the three fibres was exposed to diluted (5  $\text{L min}^{-1}$  zero air) CDS monoterpene standards for 20 min under static conditions following the previously described flushing procedure.

### 15 2.1.4 Optimisation of the sampling time

The sampling time was optimized for the 65  $\mu\text{m}$  DVB-PDMS fibre coating. Sampling times of 1, 5, 10, 20, 30, 40, 50, 60, 70, 80, 90 and 120 min were tested. Because DVB-PDMS is an adsorptive-type fibre coating, monoterpenes were extracted under non-equilibrium conditions. The sampling time was chosen as the longest practical time  
20 for which the extraction was linear such that even under the most adverse conditions (i.e. highest observed mixing ratio of target monoterpene in plant emission or in the atmosphere), the extraction does not suffer from competitive adsorption.

### 2.1.5 Calibration of monoterpenes produced by the CDS

The mixing ratios of the standards generated by diluting the CDS output with 5  $\text{L min}^{-1}$   
25 of hydrocarbon free air, were checked against the response obtained from a certified

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



calibration gas (contents listed in Table 1, National Physical Laboratory, Teddington Middlesex, UK). Relative standard deviations (RSDs) for analyses of both standards are also given in Table 1 along with the estimated theoretical mixing ratios. The difference between the measured and the theoretical concentrations range from 10 to 25% and were deemed acceptable considering the uncertainties made in the latter calculations.

### 2.1.6 Calibrations curves and limit of detection

In order to assess at which mixing ratio competitive adsorption occurs for a given fibre, mixing ratios from 3 to 455 ppbv of the diluted CDS standard gas were tested using a DVB-PDMS fibre coating and a 20 min sampling time under static conditions as previously described. A plot of the mass extracted (mass spectrometer area count) versus mixing ratio shows where competitive adsorption begins.

Diffusion rates of chemicals from vials in the CDS are a function of the vapour pressure of the pure chemical, capillary diameter, capillary length, temperature of the CDS, and the flow rates through each chamber. The ppm values achieved were dictated by the physical limits of the system as constructed and made subsequent dilution to ppt and sub-ppt levels for determination of detection limits impractical. Therefore, the detection limit of the method could not be determined simply by lowering the mixing ratios down to ppt level. For this purpose, the certified gas standard was sequentially diluted with zero-air to values as low as 2 ppt. Detection limits were measured using DVB-PDMS fiber coating and 20 min sampling time employing the static sampling procedure previously described. The detection limits of target monoterpenes are given in Table 2 and range from 4 to 20 ppt.

### 2.1.7 Storage effects

In order to evaluate the extent of adsorption of monoterpenes on the wall of the glass sampling chamber, it was filled with the monoterpene standard at 20 ppb after dilution

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



with zero air. SPME samples were then taken under static conditions at 0, 1, 2, 3, 4, 24, 48 and 72 h after the stopcocks were closed. No significant change in the MS detector response was observed, up to the third day after the chamber was filled (figure not shown). This indicates that there was no significant adsorption of monoterpenes to the walls of the glass sampling chamber.

## 2.2 Sampling from whole plant enclosure

Biogenic emissions of monoterpenes were investigated using an open, dynamic cuvette system installed inside a plant growth chamber. In other words, an air flow ran constantly through the cuvette resulting in steady state conditions mimicking natural outside conditions.

In order to validate the SPME method, monoterpene measurements were simultaneously made by a zNose-GC and also trapped on cartridges which were later analyzed by gas chromatography coupled to flame ionization detection (GC-FID). All tubes and fittings were made of PTFE and PFA Teflon. These materials are chemically inert to the relevant VOCs and feature low gas absorption and sequestration properties. Although enclosed plants can stay inside such an enclosure system for a long time without stress, physical and physiological monitoring of the plants was nevertheless performed. An overview of the components of the enclosure system is given in Fig. 2 and the components are described in detail in the following sections.

### 2.2.1 Description of cuvette chamber

Purified compressed air was used to flush the cuvette system. Compressed air was dried via passage through silica gel cartridges with a moisture indicator (orange gel, 1–3 mm, Merck, Germany). Trace gases, oil vapors, and miscellaneous particles were removed using molecular sieve beads, activated charcoal, and permanganate coated alumina (Type: IAC 630, Purafill, Germany). Glass padding and a PTFE-filter with

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion





a pore size of 5  $\mu\text{m}$  were used to remove larger particles. The silica gel and molecular sieves were periodically regenerated in a drying oven at 150 °C and 250 °C, respectively.

The cuvette was a cylindrical chamber of 0.60 m in height and 0.35 m in diameter (~58 L internal volume). In the cuvette, all surfaces in contact with sampled air are made of Teflon. The chamber consists of a skeleton of two PVC rings connected by four PVC rods, in which an FEP Teflon foil (Norton, 50  $\mu\text{m}$  thickness, Saint-Gobain Performance Plastics, Germany) is fixed. The upper side of the cuvette was closed with a lid whose surface was coated with FEP foil. The applied foil does not show interference with the trace gases tested such as monoterpenes, isoprene, organic acids or sulfur containing compounds, and exhibited only a slight absorption (5–7%) of photo-synthetic active radiation (Kesselmeier et al., 1993, 1996; Schäfer et al., 1992). The cuvette was placed inside a climate controlled chamber. In order to maintain a constant humidity in the cuvette, the air was bubbled through a vessel filled with tap water. The flow into the cuvette was controlled using a flow controller (MKS Instruments, USA) set to 5 L min<sup>-1</sup>. The air inside the cuvette was mixed by a Teflon coated propeller, driven by a magnetically coupled motor attached to the lid. The climate chamber temperature was set to 25 °C.

### 2.2.2 Micrometeorological parameters

Photosynthetically active radiation (PAR) was measured directly with a Li-Cor quantum sensor (400–700 nm; Li-Cor Biotechnology, Lincoln, NE, USA) attached to a rod of the cuvette skeleton at plant leaf level to yield realistic values for the irradiation on the leaf surface. When the light in the climate chamber was switched on, PAR values around 550  $\mu\text{mol m}^{-2} \text{s}^{-1}$  ( $\pm 25 \mu\text{mol m}^{-2} \text{s}^{-1}$ ) were obtained. Although turning off the lights resulted in an instantaneous drop in radiation, there was a short lag phase in switching on the lights. The PAR conditions were checked for every series of measurements but have been omitted in subsequent Figures for clarity in favor of white (light conditions) and grey (dark conditions) shading.

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



Humidity and air temperature was measured using a Rotronic humidity/temperature sensor (Hygrometer MP100A, Rotronic Messgeräte GmbH, Germany). The sensor was positioned in a separate chamber which was flushed with air from the cuvette using a custom-made pumping device (MPIC, Mainz, Germany), with a flow rate of 500 mL min<sup>-1</sup>. In order to avoid any radiative heating (greenhouse) effects, the chamber was enveloped with aluminum foil to reflect light.

### 2.2.3 Physiological parameters

Water vapor and CO<sub>2</sub> were analyzed online using an infrared dual channel H<sub>2</sub>O/CO<sub>2</sub> gas analyzer Li-Cor 7000 (Li-Cor Biotechnology, Lincoln, NE, USA). A custom-made pumping device (MPIC, Mainz, Germany) equipped with a built in analogue rotameter, was used to supply the analyzer with air from the cuvette with a flow rate of 500 mL min<sup>-1</sup>. Small particles were filtered from the air at the inlets of the Li-Cor using Teflon filters (Zeflour Teflon filters, 2 µm pore size, Gelman Science, USA).

### 2.2.4 Acclimation of *Quercus ilex*

The *Q. ilex* plant tested here was potted in an earth-filled bucket and positioned in the climate chamber (HPS 1500/60/S, Heraeus-Vötsch) approximately one month before the experiments started. Chamber lights were controlled to simulate day and night conditions, both lasting 12 h, with temperatures of 30 °C and 20 °C, respectively and relative humidity of 60% and 50%. The built-in lights of the chamber were switched on during this period, which irradiated the small tree with approximately 500–600 µmol m<sup>-2</sup> s<sup>-1</sup> PAR. These conditions assured reasonable, i.e. approximately natural emission of monoterpenes.

## 2.3 SPME and GC-MS analysis

The DVB-PDMS fibre coating was exposed to the air inside the whole plant enclosure by simply puncturing the Teflon foil with the SPME needle for 20 min under both light

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



and dark conditions. The experiments were performed over two separate days. The DVB-PDMS response was calibrated under dynamic conditions in the glass SPME sampling chamber (with a flow of  $5 \text{ L min}^{-1}$  out of the chamber) by diluting the CDS chemical output to 50 ppbv using  $5 \text{ L min}^{-1}$  of clean air.

Immediately after SPME sampling either from the glass SPME sampling chamber or from the cuvette, the SPME needle was introduced into the split/splitless injector of the gas chromatograph. A glass inlet liner with a narrow internal diameter (0.75 mm I.D., Supelco) was used in order to improve the GC resolution and the peak shape. Desorption was achieved in splitless mode at  $250^\circ\text{C}$  for 5 min. These settings were found to be sufficient for a quantitative desorption of all analytes studied. This was established by subjecting the analysed fibre to a second desorption and observing no carry-over peaks.

The analysis of the samples collected by SPME were conducted using the same procedure described by Yassaa and Williams (2007). Briefly, samples were separated and detected using a gas chromatograph (Agilent Technologies a GC 6890A) coupled to a mass-selective detector (Agilent Technologies MSD 5973 *inert*). The MS system was operated in electron impact mode with the following conditions: ionization potential 70 eV; source temperature  $230^\circ\text{C}$ . The MS system was operated in scan mode (30–350 u) for the identification of compounds and in SIM mode for their quantification. The enantiomeric and non-enantiomeric monoterpenes were separated using a Cyclodex-B capillary column (30 m-long, 0.256 mm I.D., 0.25  $\mu\text{m}$  film thickness) supplied by J & W Scientific (Folsom, CA, USA). The internal coating was composed of a permethylated  $\beta$ -cyclodextrin dissolved into a cyanopropyl-dimethyl polysiloxane liquid. At a constant helium carrier gas (Messer Griesheim 6.0) flow rate of  $1 \text{ ml min}^{-1}$ , the column temperature was maintained at  $40^\circ\text{C}$  for 5 min, then increased to  $160^\circ\text{C}$  at  $1.5^\circ\text{C per minute}$  as previously established by Yassaa et al. (2001). The overall analysis time was 52 min.

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

[Title Page](#)[Abstract](#)[Introduction](#)[Conclusions](#)[References](#)[Tables](#)[Figures](#)[◀](#)[▶](#)[◀](#)[▶](#)[Back](#)[Close](#)[Full Screen / Esc](#)[Printer-friendly Version](#)[Interactive Discussion](#)

## 2.4 zNose analyser

The zNose (Electronic Sensor Technology; Newbury Park, CA) is a commercially available gas chromatographic analyzer. Besides the detection of explosives and chemical warfare agents (Staples and Viswanathan, 2005), the instrument finds application in several other fields. It is used for food aroma analysis (Lammertyn et al., 2004) or for the ripeness and rot evaluation of fruits (Li et al., 2009). It has also found application for the determination of floral aroma compounds of lilac blossom (Oh et al., 2008), alarm pheromone emission from aphids (Schwartzberg et al., 2008) as well as analysis of herbivore-induced plant volatile emissions (Arimura et al., 2005; Kunert et al., 2002).

The zNose 4200 is a miniaturized high-speed gas chromatographic analyzer comprising of a pump, injector, preconcentrating trap, short GC column and a surface acoustic wave (SAW) sensor. A 6-port valve is used to trigger the trap and sensor, respectively. The system draws in a headspace sample at  $30 \text{ mL min}^{-1}$ , which is entrained onto a Tenax trap, where the compounds are adsorbed. The valve is then rotated to align the trap with the column (1 m, DB 5, ID 0.25 mm ID, 0.25  $\mu\text{m}$  film thickness). Once aligned, the trap is heated using an electrical current that vaporizes the adsorbed material. A helium carrier gas then transports the material into the column where compounds are separated. The materials exit sequentially and are channeled onto the SAW sensor.

## 2.5 VOC sampling on solid adsorbent filled cartridges and analysis with GC-FID

VOCs were sampled from the cuvette by drawing air through adsorbent filled cartridges at a constant flow rate using a custom-built, automatic cartridge sampling device, constructed at the MPIC Mainz (Kuhn et al., 2005). The system consists of a cartridge magazine comprising two sampling loops with 10 cartridge positions for sampling each, enabling sequential sampling of up to 20 cartridges. Mass flow through distinct cartridges is controlled by two miniaturized sampling valves at the inlet and at the outlet of each cartridge.

AMTD

3, 3345–3381, 2010

### Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

VOC sampling and analysis was performed in close accordance with earlier work (Kesselmeier et al., 2002). The cartridges were packed with 130 mg of Carbograph 1 (90 m<sup>2</sup> g<sup>-1</sup>) followed by 130 mg Carbograph 5 (560 m<sup>2</sup> g<sup>-1</sup>), both with a particle size of 20–40 mesh. Carbographs 1 and 5 were provided by Lara s.r.l. (Rome, Italy). The sample flow was set to 200 mL min<sup>-1</sup> and samples were taken for 15 min resulting in a sample volume of 3 L.

The analysis of the sample tubes was conducted using a GC-FID system (AutoSystem XL, Perkin-Elmer, USA). Samples were passed into the system using a thermodesorber (TurboMatrix 650, Perkin-Elmer, USA). The sensor was operated isothermally at 300 °C. A HP-1 column (100 m, 0.25 mm ID, 0.5 µm film thickness) was flushed with Helium 6.0 (Westfalen AG, Germany) at a flowrate of 1 mL min<sup>-1</sup>. The column temperature program was: 40 °C (18 min), 15 °C/min to 80 °C (40 min), 2 °C/min, 150 °C, 30 °C/min, 240 °C (5 min). The FID was calibrated using the same VOC gas standard mixture used to calibrate the CDS.

## 2.6 Application of SPME method to ambient air

The potential of SPME for monitoring VOCs in the atmosphere was tested by measuring monoterpenes present in air outside the Max Planck Institute for Chemistry (MPIC, Mainz, Germany) for two days, 14–15 August 2009. Ambient air was pumped at 5 L min<sup>-1</sup> through a 0.64 cm (1/4 inch) PTFE line from 4 m height situated directly outside the MPIC and passed through the glass SPME sampling chamber. An ozone scrubber similar to that used in the calibration experiment was installed in the front of the sampling line. After a given time, the stopcocks were closed (static mode) to capture an air sample in the chamber. The DVB-PDMS was then exposed to the enclosed air sample for 20 min. Following each air sample, the fibre was calibrated by sampling monoterpenes from the CDS and analyses were performed by GC-MS as previously described. Parallel measurements of outside air were performed by proton transfer reaction-mass spectrometry. In this technique high number densities of primary H<sub>3</sub>O<sup>+</sup> ions are continuously combined with sample air in a reaction flow drift tube maintained

### Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



at 2.2 mbar and 50 °C and ~130 Td. Any volatile organic compounds with a proton affinity higher than water react with the primary ions via proton transfer to form protonated ionic products. Primary and product ions are then sampled into a quadrupole mass spectrometer for separation and analysis. Details of the technique have been summarized most recently by Blake et al. (2009). The PTR-MS technique employs no pre-separation of chemicals. As a consequence, a variety of monoterpenes can be found with a protonated mass of  $m/z$  137 with fragments found at  $m/z$  81.

### 3 Results and discussion

#### 3.1 Comparison of SPME fibre coatings

As shown in Fig. 3 the three solid SPME fibres efficiently adsorbed the studied monoterpenes. DVB/CAR/PDMS has the highest extraction efficiency for all compounds except for myrcene and (+)-limonene. While DVB-PDMS showed a good adsorption for myrcene, (+)- $\delta$ -3-carene, (+)- $\beta$ -pinene, (–)-limonene and 1,8-cineol, CAR-PDMS showed better relative adsorption of (–)- $\alpha$ -pinene, (+)- $\alpha$ -pinene, (–)- $\beta$ -pinene and (+)-limonene. Interestingly, while (–)-limonene and (+)- $\beta$ -pinene were preferentially adsorbed on DVB-PDMS, their optical isomers, (+)-limonene and (–)- $\beta$ -pinene were better adsorbed on CAR-PDMS. This effect was even more pronounced with limonene enantiomers. The high adsorption efficiency for limonene by CAR-PDMS fibre coating was also reflected in the poor separation of limonene enantiomers (Fig. 4).

Due to the efficient adsorption of  $\alpha$ -pinene and  $\beta$ -pinene on the CAR-DVB-PDMS fibre, it becomes readily saturated with these chemicals. Competitive adsorption therefore occurred early for these species (refer to the next section for more details on competitive adsorption). This effect is most likely related to the structure and geometry of the molecules and the nature of the interactions between these molecules and the active pores of the SPME solid fibre coatings.

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

The best compromise between resolution of enantiomeric pairs and reduction of the effect of competitive adsorption was obtained using DVB-PDMS fibre coating. While DVB phase is mainly mesoporous and ideal for trapping C6–C15 analytes including monoterpenes (C10), CAR is microporous and efficiently traps C2–C6 analytes (Mani, 1999).

### 3.2 Optimisation of the sampling time

DVB-PDMS is an adsorptive SPME fibre coating and is susceptible to competitive adsorption if overloaded. Six extraction times were tested with this SPME fibre coating. The amount of monoterpene extracted is shown in Fig. 5. For the 60, 90 and 120 min extraction times, an increase in the sampling time did not lead to a significant increase in the amount of analyte extracted. This is likely due to the limited adsorptive capacity of the SPME fibre. The 1, 5, 10 and 20 min sampling times are plotted separately (Fig. 6) along with the results of a linear fit through each. For sampling times between 1 and 20 min, the linearity of the response (here expressed with respect to  $R^2$ ) ranging between 0.96 and 0.99, was deemed acceptable for use in quantification. Within the 20 min sampling time, there was no apparent displacement of the compounds and 20 min was therefore chosen as the optimum sampling time for DVB-PDMS fibre coating to quantify monoterpenes.

As can be seen from Fig. 5,  $\alpha$ -pinene,  $\beta$ -pinene and 1,8-cineol reach equilibrium earlier than other monoterpene species due to the displacement effect. Such competition on DVB-PDMS can be avoided by limiting sampling to 20 min in situations where the mixing ratio of any given monoterpene does not exceed 100 ppb. It should be noted that in using a DVB-PDMS SPME fibre under non-equilibrium conditions, the total response represents the time averaged mixing ratio encountered by the fibre over the course of the sampling period.

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion





### 3.3 Calibration curves and detection limits

The quantification of target monoterpenes was based on the calibration curves for the optimized 20 min sampling time. The mixing ratio ranges of the monoterpenes are shown in Table 3. Except  $\alpha$ -pinene, which responds linearly only up to around 53 ppb, most chemicals respond linearly up to 100 ppb level. Limonene responded linearly up to 300 ppb. These results showed how  $\alpha$ -pinene is affected by the competitive adsorption at high mixing ratios, a fact which should be considered during calibration. Since the maximum measured mixing ratios of the most abundant monoterpene are around a few ppb in the atmosphere and rarely exceed 50 ppb in plant sampling (Kesselmeier and Staudt, 1999), the limit of linearity range measured in this study is quite acceptable.

The limit of detection showed in Table 2 ranged from 4 ppt for  $\alpha$ -pinene to 20 ppt for *cis*- $\beta$ -ocimene. The range of detection limits is quite reasonable and comparable to those obtained in sampling with solid sorbent and thermal desorption.

### 3.4 Application to whole plant enclosure

A *Q. ilex* plant was investigated using SPME, zNose and cartridge VOC sampler side by side. Figure 7 shows the results of the zNose measurements. The plotted values are the signal peak area of measurements with sample times of 2 min at 40 °C sensor temperature. Four signals could be observed, which showed an obvious light dependency with Kovats Indices (KIs) of 988, 1031, 1207 and 1612. Shortly after the light was switched on, emissions of the four compounds could be observed and reached maximum values after approximately 45 min. When the light was turned off, the emissions dropped below the detection limit within approximately 15 min. Using the observed KIs, signals can be tentatively assigned to limonene (KI of 1031) and  $\beta$ -myrcene (KI of 988).

On the other hand, Fig. 8 shows an example of the reconstructed mass chromatographic profile of monoterpenes obtained by analyzing the SPME sample collected from the gaseous emission of *Q. ilex* under light and dark (overlaid in blue) condition on a  $\beta$ -cyclodextrin column using single ion mode (SIM). The chromatographic profile

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion





depicted in Fig. 8 was generated by plotting the current produced by ions with  $m/z$  93, a fragment quite specific to monoterpene hydrocarbons that, in many cases, is the base peak of the mass spectrum. The chromatogram shows only two peaks that are attributed to  $\beta$ -myrcene and (–)-limonene. By comparing the ratios of  $\beta$ -myrcene to limonene between SPME GC-MS and zNose, we conclude that the peaks KI 1031 and KI 988 obtained from zNose during these experiments are very likely limonene and  $\beta$ -myrcene, respectively.

The SPME (Table 4) and VOC analysis by GC-FID (Fig. 9) showed  $\beta$ -myrcene as the most prominent chemical emitted, followed by limonene. These two chemicals made up more than 95% of the total monoterpene emissions. Very small fractions of  $\alpha$ -pinene,  $\beta$ -pinene and sabinene, as well as a few other monoterpenes, made up the residual percentages.

While the derived gas exchange rates for limonene from SPME revealed average values of  $220 \text{ ng m}^{-2} \text{ s}^{-1}$ , the ones from cartridge GC-FID were around  $240 \text{ ng m}^{-2} \text{ s}^{-1}$  and those from zNose ranged from 250 to 300 (see, Fig. 10). The gas exchange rates from  $\beta$ -myrcene were  $590 \text{ ng m}^{-2} \text{ s}^{-1}$  derived from SPME and  $610 \text{ ng m}^{-2} \text{ s}^{-1}$  obtained from cartridge GC-FID. However, the comparison could not be made with the zNose, because no response factor for myrcene was available.

Previous studies have shown  $\alpha$ -pinene and  $\beta$ -pinene as the most prominent compounds emitted by *Q. ilex* with 9,9% limonene and 6,7%  $\beta$ -myrcene found in lesser amounts (Kesselmeier et al., 1996; Yassaa et al., 2001). Three main emission types of *Q. ilex* have been distinguished in literature. Of investigated populations, 70.5% emitted mainly  $\alpha$ -pinene and  $\beta$ -pinene (60–80%), 21.2% produced higher proportions of limonene (30–80%), and 8.2% had higher proportions of  $\beta$ -myrcene (60–90%) in their emissions (Staudt et al., 2001). The observed monoterpene emissions were therefore not exceptional and simply reflect the different emitter types of this species.

The exchange rates of the potted individual enclosed in the cuvette showed similar results for  $\beta$ -myrcene by the three techniques and by SPME and GC-FID for limonene. The observed emission rates were comparable and fit into the range of reported to-

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

tal monoterpene emission rates of *Q. ilex*. This value is generally between 10 and 1600 ng m<sup>-2</sup> s<sup>-1</sup> depending on the individual plant, season, and measurement conditions (Staudt et al., 2001).

It should be noted that the studied *Q. ilex* individual showed a strong enantiomeric preference (almost 100% (–)-limonene) compared to the previously investigated species by Yassaa et al. (2001) where (–)-limonene and (+)-limonene accounted 41 and 59%, respectively.

### 3.5 Application to ambient air

The mixing ratios of monoterpenes in the ambient air of outside of MPIC (Mainz, Germany) at 4 m height were determined over two days. The samples were taken by filling the 5 L glass SPME sampling chamber using a pump. Both ambient air and calibration gas were sampled with SPME under static conditions with the two valves closed and a sample time of 20 min. Parallel measurements were performed with proton transfer reaction-mass spectrometry with *m/z* 137 used to quantify monoterpenes. Since PTR-MS is not selective for monoterpenes, only the mixing ratios of total monoterpenes could be determined for comparison.

Table 5 reports the min, max and average of 2 days (14–15 August 2009) of measurements of monoterpenes in ambient air outside of MPIC (Mainz, Germany). The total mixing ratios of monoterpenes by PTR-MS, 600 pptv was close to that determined by SPME (425 pptv as the sum of the identified and quantified monoterpenes). This provides additional support for the applicability of the SPME system to the measurement of monoterpenes in the atmosphere.

## 4 Conclusions

Using optimised 20 min sampling times and conditions with a DVB-PDMS fibre coating, SPME has shown good potential for studying BVOCs both in a plant chamber and in

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



the atmosphere. Emission rates of monoterpenes from a potted individual of *Q. Ilex* determined by SPME are comparable to rates found by the zNose and through use of solid sorbent and GC-FID.

Although it has been shown that the zNose can qualitatively measure the emission potential and kinetics of plant volatile emission, chemical speciation must be confirmed using alternative methods. SPME represents the technique of choice for this purpose. Compared to cartridge based adsorption and thermal desorption techniques, SPME is cheap, simple, fast, sensitive. It is highly mobile, completely eliminates the use of organic solvents, and does not require the use of an expendable cryogen in the analysis. This easily automated method when combined with a portable and fast GC. These qualities should provide excellent opportunities for studying emissions of organic volatiles from a variety of different plant species, both in plant chambers and in ambient air.

*Acknowledgements.* The authors would like to thank Uwe Kuhn, Tamara Dindorf, Thomas Kenntner, Frank Helleis and Nina-Maria Knothe from Biogeochemistry Department, Max-Planck Institute for Chemistry for their technical support in constructing the permeation device.

The service charges for this open access publication have been covered by the Max Planck Society.

## References

- Arimura, G., Kost, C., and Boland, W.: Herbivore-induced, indirect plant defences, *Biochim. Biophys. Acta*, 1734, 91–111, 2005.
- Baldwin, I. T., Halitschke, R., Paschold, A., von Dahl, C. C., and Preston, C. A.: Volatile signaling in plant-plant interactions: “talking trees” in the genomics era, *Science*, 311, 812–815, 2006.
- Baker, B. and Sinnott, M.: Analysis of sesquiterpene emissions by plants using solid phase microextraction, *J. Chromatogr. A*, 1216, 8442–8451, 2009.
- Bicchi, C., Drigo, S., and Rubiolo, P.: Influence of fibre coating in the headspace solid-phase microextraction-gas chromatographic analysis of aromatic and medicinal plants, *J. Chromatogr. A*, 892, 469–485, 2000.

## Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

## Abstract

## Introduction

## Conclusions

## References

## Tables

## Figures



▶

[Back](#)

Close

Full Screen / Esc

Printer-friendly Version

## Interactive Discussion



- Blake, R. S., Monks, P. S., and Ellis, A. M.: Proton-transfer reaction mass spectrometry, *Chem. Rev.*, 109, 861–896, doi:10.1021/cr800364q, 2009.
- Croteau, R.: Biosynthesis and catabolism of monoterpenoids, *Chem. Rev.*, 87, 929–954, 1987.
- Guenther, A., Hewitt, C. N., Erickson, D., Fall, R., Geron, C., Graedel, T., Harley, P., Klinger, L.,  
5 Lerda, M., McKay, W. A., Pierce, T., Scholes, B., Steinbrecher, R., Tallamraju, R., Taylor, J.,  
and Zimmerman, P.: A global-model of natural volatile organic compound emissions, *J. Geo-  
phys. Res.*, 100, 8873–8892, 1995.
- Helmig, D., Reevermann, T., Pollmann, J., Kaltschmidt, O., Hernandez, A. J., Bocquet, F., and  
David, D.: Calibration system and analytical considerations for quantitative sesquiterpene  
10 measurements in air, *J. Chromatogr. A*, 1002, 193–211, 2003.
- Kesselmeier, J., Meixner, F. X., Hofmann, U., Ajavon, A.-L., Leimbach, S., and Andreae, M. O.:  
Reduced sulfur compound exchange between the atmosphere and tropical tree species in  
Southern Cameroon, *Biogeochemistry*, 23, 23–45, 1993.
- Kesselmeier, J., Schäfer, L., Ciccioli, P., Brancaleoni, E., Cecinato, A., Frattoni, M., Foster, P.,  
15 Jacob, V., Denis, J., Fugit, J. L., Dutaur, L., and Torres, L.: Emission of monoterpenes and  
isoprene from a Mediterranean oak species *Quercus ilex* L. measured within the BEMA  
(Biogenic Emissions in the Mediterranean Area) project, *Atmos. Environ.*, 30, 1841–1850,  
1996.
- Kesselmeier, J. and Staudt, M.: Biogenic volatile organic compounds (VOC): an overview on  
20 emission, physiology and ecology, *J. Atmos. Chem.*, 33, 23–88, 1999.
- Kesselmeier, J., Kuhn, U., Rottenberger, S., Biesenthal, T., Wolf, A., Schebeske, G., An-  
dreae, M. O., Ciccioli, P., Brancaleoni, E., Frattoni, M., Oliva, S. T., Botelho, M. L.,  
Silva, C. M. A., and Tavares, T. M.: Concentrations and species composition of atmospheric  
volatile organic compounds (VOC) as observed during the wet and dry season in Rondônia  
25 (Amazonia), *J. Geophys. Res.*, 107, 8053, 2002.
- Kuhn, U., Dindorf, T., Ammann, C., Rottenberger, S., Guyon, P., Holzinger, R., Ausma, S.,  
Kenntner, T., Helleis, F., and Kesselmeier, J.: Design and field application of an automated  
cartridge sampler for VOC concentration and flux measurements, *J. Environ. Monitor*, 7,  
568–76, 2005.
- Kunert, M., Biedermann, A., Koch, T., and Boland, W.: Ultrafast sampling and analysis of plant  
30 volatiles by a hand-held miniaturised GC with pre-concentration unit: kinetic and quantitative  
aspects of plant volatile production, *J. Sep. Sci.*, 25, 677–684, 2002.
- Lammertyn, J., Veraverbeke, E. A., and Irudayaraj, J.: zNose™ technology for the classifica-

- tion of honey based on rapid aroma profiling, *Sensor. Actuat. B-Chem.*, 98, 54–62, 2004.
- Li, Z., Wang, N., Raghavan, G. V., and Vigneault, C.: Ripeness and rot evaluation of “Tommy Atkins” mango fruit through volatiles detection, *J. Food Eng.*, 91, 319–324, 2009.
- Lough, W. J. and Wainer, I. W.: *Chirality in Natural and Applied Science*, CRC Press, Blackwell Publishing, Oxford, 2002.
- Mani, V.: Properties of commercial SPME coatings, in: *Applications of Solid Phase Microextraction*, edited by: Pawliszyn, J., The Royal Society of chemistry, Cambridge, UK, 57–108, 1999.
- Müller, J. F.: Geographical distribution and seasonal variation of surface emissions and deposition velocities of atmospheric trace gases, *J. Geophys. Res.*, 97, 3787–3804, 1992.
- Oh, S. Y., Ko, J. W., Jeong, S.-Y., and Hong, J.: Application and exploration of fast gas chromatography-surface acoustic wave sensor to the analysis of thymus species, *J. Chromatogr. A*, 1205, 117–127, 2008.
- Pollmann, J., Ortega, J., and Helmig, D.: Analysis of atmospheric sesquiterpenes: sampling losses and mitigation of ozone interference, *Environ. Sci. Technol.*, 39, 9620–9629, 2005.
- Schäfer, L., Kesselmeier, J., and Helas, G.: Formic and acetic acid emission from conifers measured with a “cuvette” technic, in: *CeC Air Pollution Research 39: Field Measurements and Interpretation of Species Related to Photooxidants and Acid Deposition*, edited by: Angeletti, G., Beilke, S., and Slanina, J., Eur. Comm., Brussels, 319–323, 1992.
- Schwartzberg, E., Kunert, G., Stephan, C., David, A., Röse, U., Gershenzon, J., Boland, W., and Weisser, W.: Real-time analysis of alarm pheromone emission by the pea aphid (*Acyrtosiphon pisum*) under predation, *J. Chem. Ecol.*, 34, 76–81, 2008.
- Staples, E. and Viswanathan, S.: Ultrahigh-speed chromatography and virtual chemical sensors for detecting explosives and chemical warfare agents, *IEEE Sens. J.*, 5, 622–631, 2005.
- Staudt, M., Mandl, N., Joffre, R., and Rambal, S.: Intraspecific variability of monoterpene composition emitted by *Quercus ilex* leaves, *Can. J. Forest Res.*, 31, 174–180, 2001.
- Stokes, G. Y., Chen, E. H., Buchbinder, A. M., Paxton, W. F., Keeley, A., and Geiger, F. M.: Atmospheric heterogeneous stereochemistry, *J. Am. Chem. Soc.*, 131, 13733–13737, 2009.
- Williams, J.: Organic trace gases: an overview, *Environ. Chem.*, 1, 125–136, 2004.
- Williams, J., Yassaa, N., Bartenbach, S., and Lelieveld, J.: Mirror image hydrocarbons from Tropical and Boreal forests, *Atmos. Chem. Phys.*, 7, 973–980, doi:10.5194/acp-7-973-2007, 2007.
- Unsicker, S. B., Kunert, G., and Gershenzon, J.: Protective perfumes: the role of vegetative

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

I◀

▶I

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



volatiles in plant defense against herbivores, Curr. Opin. Plant Biol., 12, 479–485, 2009.

Yassaa, N., Meklati, B. Y., and Cecinato, A.: Evaluation of monoterpenic biogenic volatile organic compounds in ambient air around *Eucalyptus globulus*, *Pinus halepensis* and *Cedrus atlantica* trees growing in Algiers city area by chiral and achiral capillary gas chromatography, Atmos. Environ., 34, 2809–2816, 2000.

Yassaa, N., Brancaleoni, E., Frattoni, M., and Cicciooli, P.: Trace level determination of enantiomeric monoterpenes in terrestrial plant emission and in the atmosphere using a  $\beta$ -cyclodextrin capillary column coupled with thermal desorption and mass spectrometry, J. Chromatogr. A, 915, 185–197, 2001.

Yassaa, N. and Williams, J.: Analysis of enantiomeric and non-enantiomeric monoterpenes in plant emissions using portable dynamic air sampling/solid-phase microextraction (PDAS-SPME) and chiral gas chromatography/mass spectrometry, Atmos. Environ., 39, 4875–4884, 2005.

Yassaa, N. and Williams, J.: Enantiomeric monoterpene emissions from natural and damaged Scots pine in a boreal coniferous forest measured using solid-phase microextraction and gas chromatography/mass spectrometry, J. Chromatogr. A, 1141, 138–144, 2007.

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

**Table 1.** Mixing ratios of monoterpenes from the standard gas bottle, measured from CDS and calculated from the weight loss of CDS tubes.

	Gas standard bottle	CDS			
	Mixing ratios (ppb)	Channel	Flow rate (SCCM)	Measured mixing ratios (ppm) for CDS	Theoretical mixing ratios (ppm) from weight loss of CDS tubes
(1 <i>S</i> )-(-)- $\alpha$ -pinene	2.00±0.06	1	58.35	1.75±0.07	2.06
(1 <i>S</i> )-(-)- $\beta$ -pinene	1.98±0.06	1		2.01±0.07	2.51
(4 <i>S</i> )-(-)-limonene	2.05±0.06	1		1.56±0.08	1.75
$\beta$ -myrcene	1.92±0.10	2	53.64	2.27±0.11	1.66
(+)- $\delta$ -3-carene	1.99±0.06	2		1.24±0.07	1.62
1,8-cineol	2.18±0.11	2		1.49±0.12	1.65
(1 <i>R</i> )-(+)- $\alpha$ -pinene	2.00±0.06	3	54.64	1.84±0.07	2.17
(1 <i>R</i> )-(+)- $\beta$ -pinene	1.98±0.06	3		2.14±0.07	1.87
(4 <i>R</i> )-(+)-limonene	2.05±0.06	3		2.10±0.08	1.81

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



**Table 2.** Detection limits of monoterpenes sampled from glass sampling chamber using DVB-PDMS fibre coating and gas standard bottle.

Chemical	Mixing ratios (ppt)
(1 <i>S</i> )-(-)- $\alpha$ -pinene	4
(1 <i>R</i> )-(+)- $\alpha$ -pinene	4
$\beta$ -myrcene	5
(1 <i>S</i> ,5 <i>S</i> )-(-)-sabinene	10
$\delta$ -3-carene	6
(1 <i>R</i> )-(+)- $\beta$ -pinene	5
(1 <i>S</i> )-(-)- $\beta$ -Pinene	5
Cis- $\beta$ -ocimene	20
<i>p</i> -cymene	4
(-)-limonene	5
1,8-cineol	8



**Table 3.** Calibration curves with correlation coefficients of monoterpenes sampled from the glass sampling chamber using DVB-PDMS and CDS at different mixing ratios.

Mixing ratio range (ppb)	3–15	3–21	3–35	3–53	3–105	3–157
(+)- $\alpha$ -pinene	$y=188\,800x+308\,296$ $R^2=0.9992$	$y=168\,820x+430\,482$ $R^2=0.9896$	$y=128\,381x+808\,543$ $R^2=0.9505$	$y=90\,543x+1E+06$ $R^2=0.876$	$y=47\,522x+2E+06$ $R^2=0.7348$	$y=30\,533x+3E+06$ $R^2=0.6692$
Mixing ratio range (ppb)	12–60	12–80	12–134	12–201	12–402	12–602
$\beta$ -Myrcene	$y=86\,728x-260\,317$ $R^2=0.9989$	$y=85\,038x-220\,729$ $R^2=0.9992$	$y=75\,910x+106\,226$ $R^2=0.993$	$y=60\,799x+839\,360$ $R^2=0.9611$	$y=38\,778x+2E+06$ $R^2=0.8941$	$y=28\,249x+3E+06$ $R^2=0.8618$
Mixing ratio range (ppb)	6–33	6–44	6–73	6–110	6–219	6–329
(+)- $\delta$ -3-carene	$y=133\,670x-28\,862$ $R^2=0.9994$	$y=123\,944x+95\,667$ $R^2=0.995$	$y=105\,858x+449\,699$ $R^2=0.9846$	$y=81\,814x+1E+06$ $R^2=0.9429$	$y=50\,130x+2E+06$ $R^2=0.8654$	$y=34\,250x+3E+06$ $R^2=0.804$
Mixing ratio range (ppb)	4–22	4–29	4–49	4–73	4–147	1–220
(+)- $\beta$ -pinene	$y=14\,3031x+123\,372$ $R^2=1$	$y=12\,8183x+250\,549$ $R^2=0.9904$	$y=10\,3050x+579\,660$ $R^2=0.9685$	$y=75\,813x+1E+06$ $R^2=0.9106$	$y=43\,009x+2E+06$ $R^2=0.8048$	$y=28\,215x+2E+06$ $R^2=0.7344$
Mixing ratio range (ppb)	9–45	9–61	9–101	9–152	9–303	9–455
(+)-limonene	$y=146\,480x-201\,076$ $R^2=0.9989$	$y=139\,133x-71\,141$ $R^2=0.9975$	$y=122\,308x+38\,3731$ $R^2=0.9902$	$y=97\,167x+1E+06$ $R^2=0.9568$	$y=62\,103x+3E+06$ $R^2=0.8929$	$y=44\,183x+5E+06$ $R^2=0.8474$
Mixing ratio range (ppb)	7–35	7–50	7–83	7–125	7–250	7–375
1,8-cineol	$y=74\,551x+44\,304$ $R^2=0.9997$	$y=69\,743x+114\,321$ $R^2=0.9964$	$y=59\,532x+341\,675$ $R^2=0.9849$	$y=46\,339x+739\,537$ $R^2=0.946$	$y=28\,534x+2E+06$ $R^2=0.8696$	$y=19\,959x+2E+06$ $R^2=0.8219$

**Quantitative analysis  
of chiral  
monoterpenes using  
SPME**

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion



**Table 4.** Average gas exchange rates ( $\text{ng m}^{-2} \text{s}^{-1}$ ) of monoterpenes from the emission of potted individual of *Q. ilex* obtained with SPME with DVB-PDMS fibre at 20 min sampling during light and dark conditions compared to those recorded by solid sorbent GC-FID and zNose.

	SPME Light off	SPME Light on	Solid sorbent GC-FID	zNose
(1 <i>S</i> )-(-)- $\alpha$ -pinene	u.d.l.*	1.5		
(1 <i>R</i> )-(+)- $\alpha$ -pinene	u.d.l.	0.6		
$\beta$ -myrcene	u.d.l.	592	610	n.q.**
(1 <i>S</i> ,5 <i>S</i> )-(-)-sabinene	u.d.l.	5.5		
(1 <i>S</i> )-(-)-camphene	u.d.l.	0.2		
(1 <i>R</i> )-(+)- $\beta$ -pinene	u.d.l.	0.5		
(1 <i>S</i> )-(-)- $\beta$ -pinene	u.d.l.	0.5		
Cis- $\beta$ -ocimene	u.d.l.	0.5		
(4 <i>S</i> )-(-)-limonene	0.6	223	240	250–300
(4 <i>R</i> )-(+)-limonene	u.d.l.	1.0		
1,8-cineol	u.d.l.	0.3		

\* u.d.l.: under detection limit.

\*\* n.q.: detected but not quantified.

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

**Table 5.** Average monoterpene mixing ratios (pptv) of monoterpenes present in the ambient air outside MPIC (Mainz, Germany) sampled with SPME, 14–15 August 2009.

	Min	Max	Average $\pm$ RSD%
(1 <i>S</i> )-(-)- $\alpha$ -pinene	60	65	64 $\pm$ 5
(1 <i>R</i> )-(+)- $\alpha$ -pinene	49	56	53 $\pm$ 6
$\beta$ -myrcene	22	24	24 $\pm$ 3
(+)- $\delta$ -3-carene	45	51	49 $\pm$ 6
(1 <i>R</i> )-(+)- $\beta$ -pinene	32	45	40 $\pm$ 18
(1 <i>S</i> )-(-)- $\beta$ -pinene	52	62	58 $\pm$ 9
(4 <i>S</i> )-(-)-limonene	31	36	35 $\pm$ 8
(4 <i>R</i> )-(+)-limonene	36	40	39 $\pm$ 6
1,8-cineol	57	61	59 $\pm$ 3

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶


Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

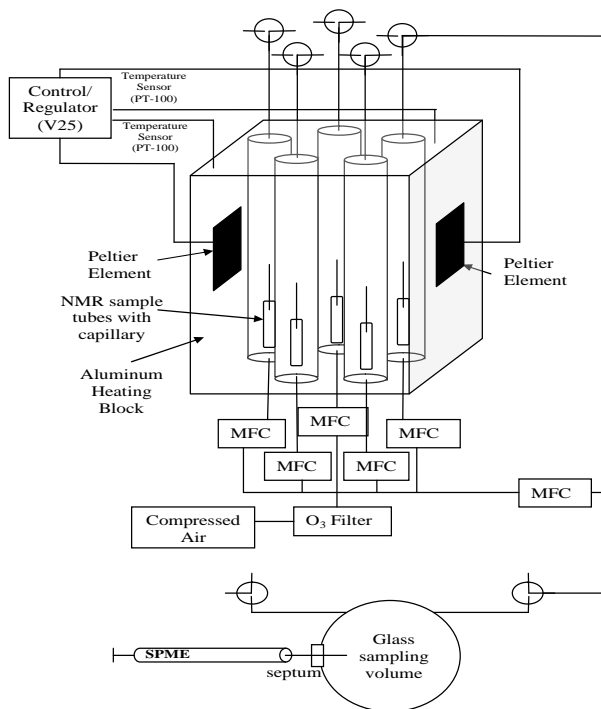

**Fig. 1.** Schematic representation of capillary diffusion system (CDS) and glass sampling chamber. Aluminum block is maintained at 27 °C by either heating or cooling the peltier elements. Entire aluminum block is surrounded by a thick layer of insulation material and peltiers are connected to large cooling vanes to effectively dissipate heat. Tubing leading to the glass chambers is pre-heated to equilibrate air temperature prior to entry into the glass chamber. MFC: mass flow controller. V25: a MPI custom built control/regulator.

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close


Full Screen / Esc

Printer-friendly Version

Interactive Discussion

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

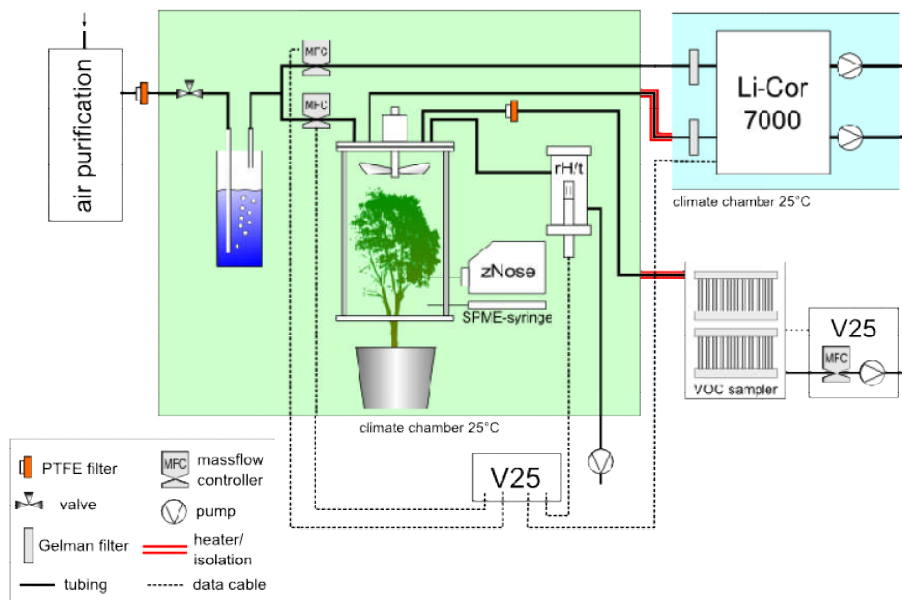

**Fig. 2.** Overview of the cuvette system. The SPME and zNose are attached to the cuvette system by simply puncturing the Teflon foil with their respective needles.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close


Full Screen / Esc

Printer-friendly Version

Interactive Discussion


# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

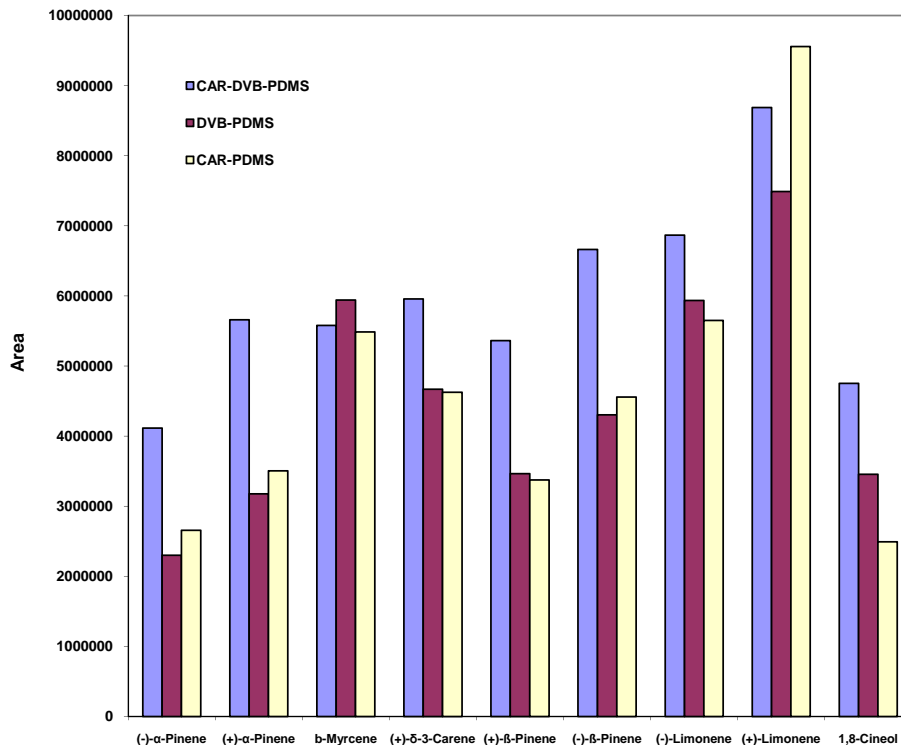

**Fig. 3.** A comparison of the extraction efficiencies of DVB-PDMS, CAR-PDMS and CAR-DVB-PDMS for monoterpenes.


[Title Page](#)
[Abstract](#)
[Introduction](#)
[Conclusions](#)
[References](#)
[Tables](#)
[Figures](#)
[◀](#)
[▶](#)
[◀](#)
[▶](#)
[Back](#)
[Close](#)
[Full Screen / Esc](#)
[Printer-friendly Version](#)
[Interactive Discussion](#)

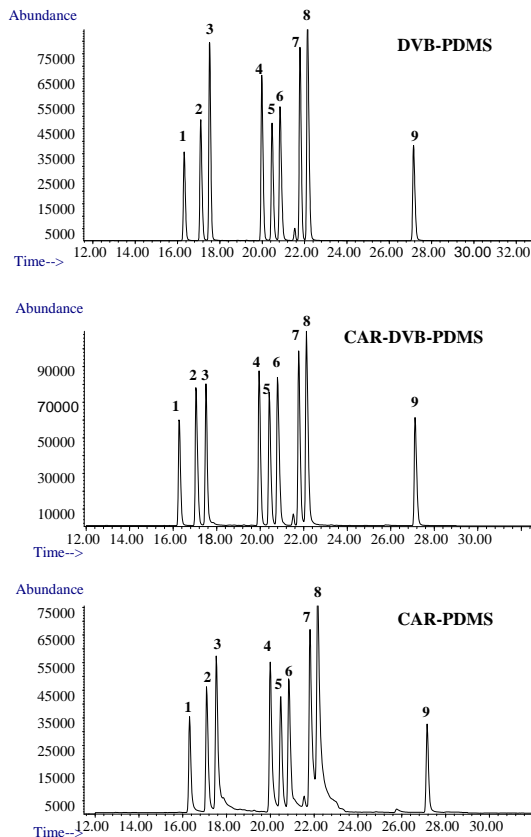

**Fig. 4.** Reconstructed mass chromatographic profiles ( $m/z$  93) of monoterpenes sampled with three different SPME fibre coatings with 20 min sampling time. (1):  $(-)\text{-}\alpha$ -pinene, (2):  $(+)\text{-}\alpha$ -pinene, (3):  $\beta$ -myrcene, (4):  $(+)\text{-}\delta$ -3-carene, (5):  $(+)\text{-}\beta$ -pinene, (6):  $(-)\text{-}\beta$ -pinene, (7):  $(-)$ -limonene, (8):  $(+)\text{-}$ limonene and (9): 1,8-cineol.

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

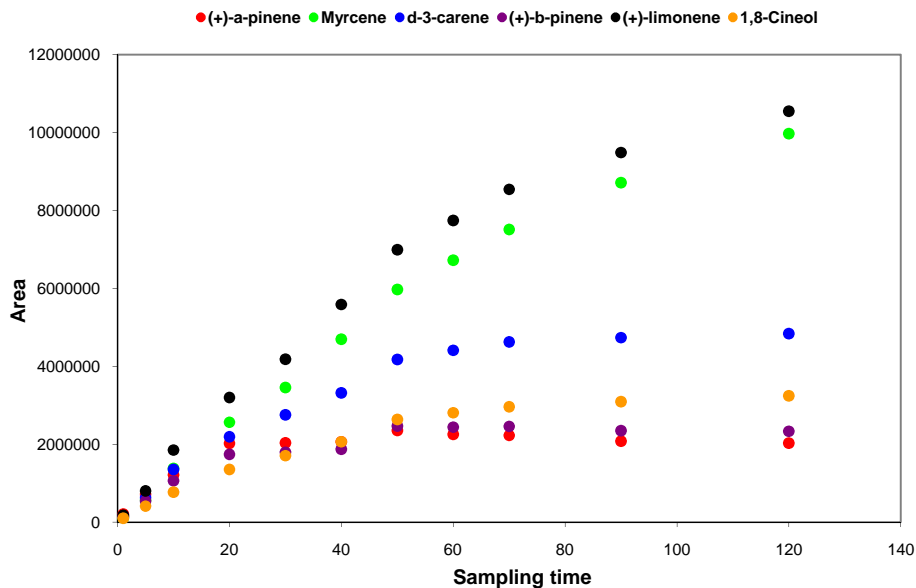


**Fig. 5.** Extraction efficiencies of monoterpenes in 1, 5, 10, 20, 30, 40, 50, 60, 70, 80, 90 and 120 min sampling times using DVB-PDMS.

[Title Page](#)
[Abstract](#)
[Introduction](#)
[Conclusions](#)
[References](#)
[Tables](#)
[Figures](#)
[◀](#)
[▶](#)
[◀](#)
[▶](#)
[Back](#)
[Close](#)
[Full Screen / Esc](#)
[Printer-friendly Version](#)
[Interactive Discussion](#)



# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

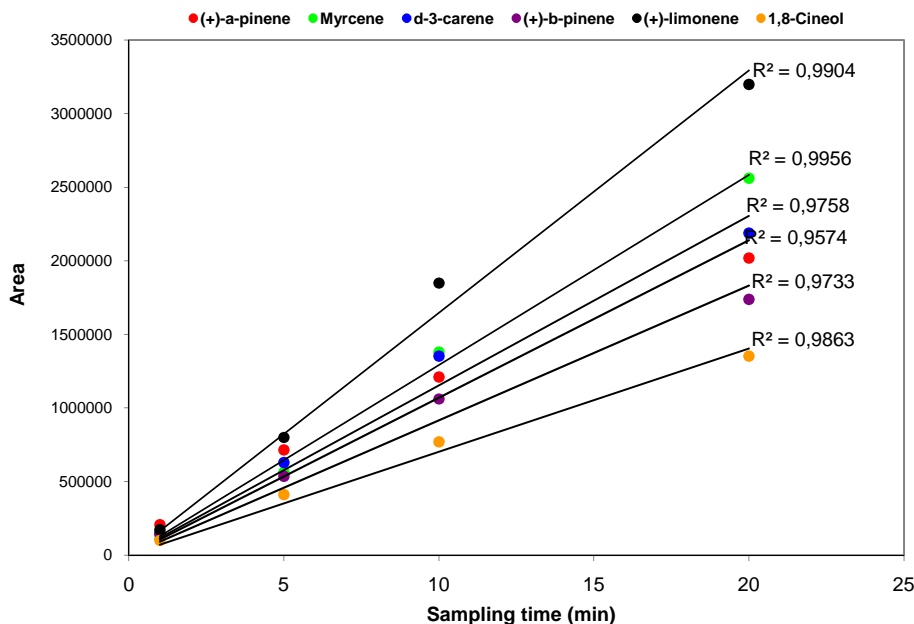


**Fig. 6.** Extraction efficiencies of monoterpenes in 1, 5, 10 and 20 min sampling times and correlation coefficients using DVB-PDMS.

[Title Page](#)
[Abstract](#)
[Introduction](#)
[Conclusions](#)
[References](#)
[Tables](#)
[Figures](#)
[◀](#)
[▶](#)
[◀](#)
[▶](#)
[Back](#)
[Close](#)
[Full Screen / Esc](#)
[Printer-friendly Version](#)
[Interactive Discussion](#)


# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

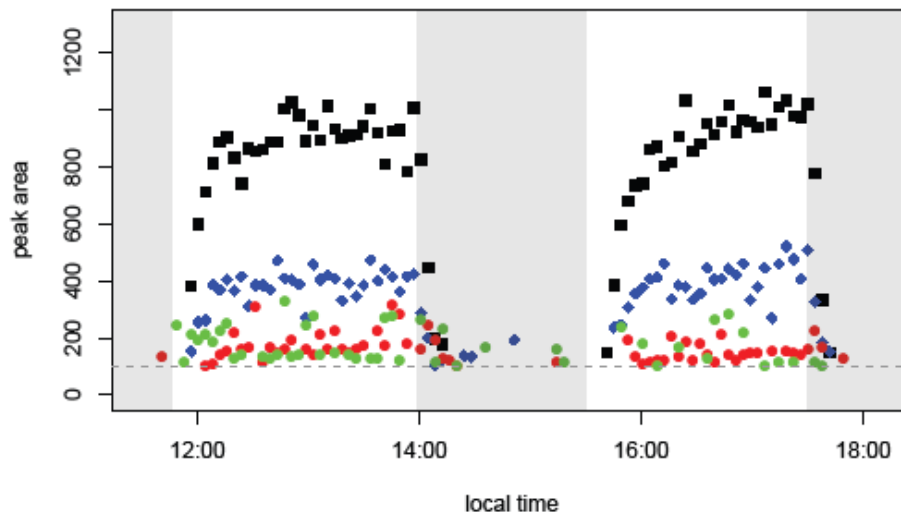

**Fig. 7.** Emissions of the potted *Q. ilex*, measured in the cuvette with sample times of 3 min at 40 °C sensor temperature. The dashed line indicates the peak detection limit. i KI 988, m KI 1031, Y KI 1207, Y KI 1612.

[Title Page](#)
[Abstract](#)
[Introduction](#)
[Conclusions](#)
[References](#)
[Tables](#)
[Figures](#)
[◀](#)
[▶](#)
[◀](#)
[▶](#)
[Back](#)
[Close](#)
[Full Screen / Esc](#)
[Printer-friendly Version](#)
[Interactive Discussion](#)

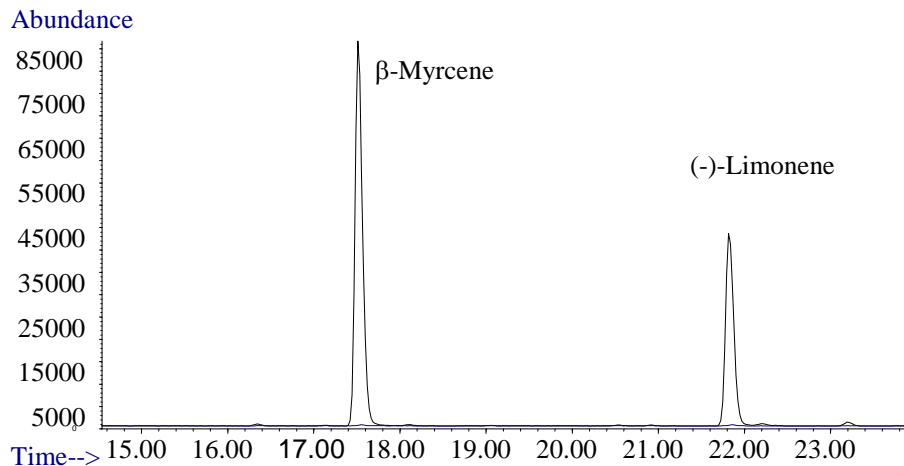

**Fig. 8.** Reconstructed mass chromatographic profiles ( $m/z$  93) of monoterpenes present in the emission of potted *Q. ilex* during light and dark (overlaid in blue) conditions.

**Quantitative analysis  
of chiral  
monoterpenes using  
SPME**

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶

Back

Close


Full Screen / Esc

Printer-friendly Version

Interactive Discussion

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

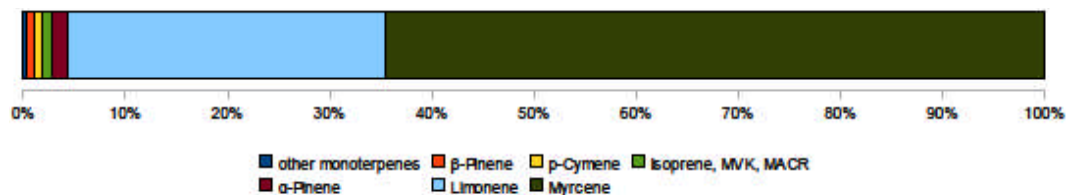

**Fig. 9.** Composition of monoterpenes emitted by the potted individual of *Q. ilex* according to solid sorbent GC-FID measurements.

[Title Page](#)
[Abstract](#)
[Introduction](#)
[Conclusions](#)
[References](#)
[Tables](#)
[Figures](#)
[◀](#)
[▶](#)
[◀](#)
[▶](#)
[Back](#)
[Close](#)
[Full Screen / Esc](#)
[Printer-friendly Version](#)
[Interactive Discussion](#)

# Quantitative analysis of chiral monoterpenes using SPME

N. Yassaa et al.

Title Page

Abstract

Introduction

Conclusions

References

Tables

Figures

◀

▶

◀

▶


Back

Close

Full Screen / Esc

Printer-friendly Version

Interactive Discussion

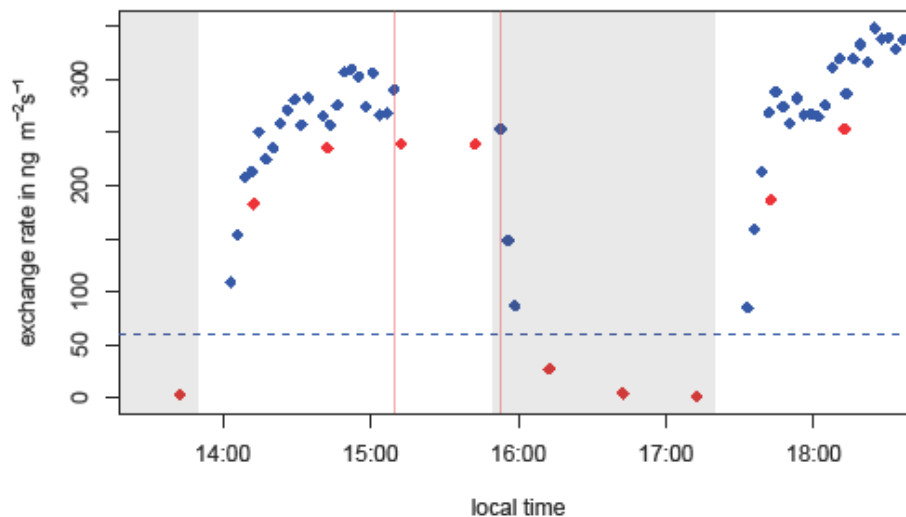

**Fig. 10.** Limonene calculated gas exchange rates for the zNose data (m), KI 1031 compared to rates obtained by solid sorbent GC-FID data (m). The dashed line indicates the peak detection limit of zNose.